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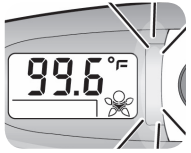
49545 & 49657 User Guide Exchangeable Tip 3-in-1 Thermometer


! WARNING:
BEFORE USING DIGITAL THERMOMETER, PLEASE READ CAREFULLY.
• This is **NOT** an ear thermometer. Take temperature measurements Underarm, Orally or Rectally.
• **NEVER** allow a child to play with these products. They are not toys. They are sensitive electronic devices and may be damaged if not used properly under parental supervision.
• **ALWAYS** supervise your child when using thermometer. **NEVER** allow your child to use it alone or without close parental supervision.
• Avoid sudden impact and dust.

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Know Your Unit

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
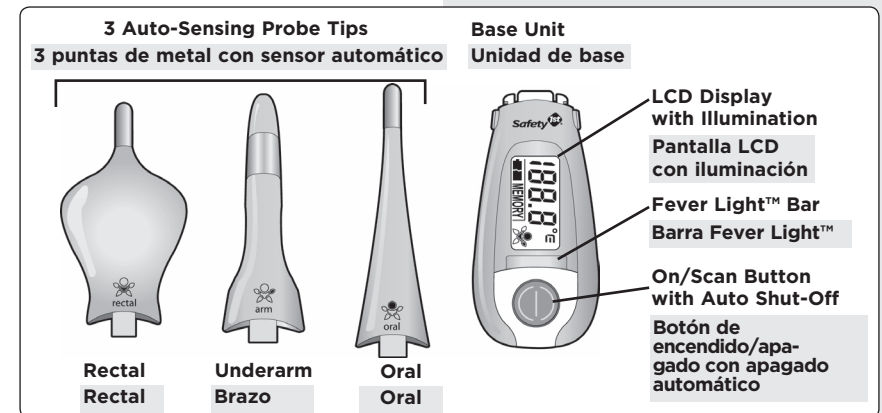
To Use Tips

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- To attach a tip, select oral, rectal, or underarm tip. Align icon to the front (above the Safety 1st logo) and plug tip into base unit (**Figure 1**). Tip fits in one direction only. The mode is automatically sensed and will display accordingly (**Figure 1**).
- To remove tip, press buttons on both sides and pull tip off (**Figure 2**).

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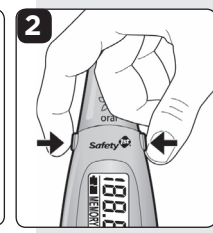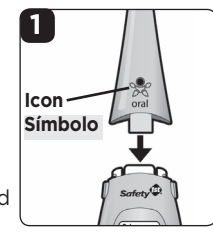Important Information

Información importante

- Read these instructions carefully to ensure accurate temperature readings. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.
- High, prolonged fever requires medical attention especially for young children. Please contact your physician.
- Do not allow children to take their temperatures unattended. This thermometer is not a toy and should be kept out of a child's reach when not in use.
- Do not allow children to walk or run during temperature taking. The probe could present a hazard.
- Oral, underarm, and rectal temperatures vary by approximately 1.0°F (0.5°C).
- Do not return this product to the place of purchase. If any parts are missing, email consumer@djgusa.com, call Consumer Relations at (800) 544-1108, or fax at (800) 207-8182. You can also visit our website at www.safety1st.com. Have the model number (49545 or 49657) and date code (manufacture date located on back of thermometer) ready.

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Temperature Measurements Range

Many people believe that normal body temperature is 98.6°F (37.0°C). However, medical research has shown that temperature varies from person to person. Factors affecting body temperature include age, gender, physical activity, environmental temperature, clothing, medications, and time of day. Oral, rectal, and underarm temperatures vary by approximately 1.0°F/0.5°C.

Method	Average Normal Temperature Range	Fever*
Oral	97.0°F to 99.6°F (36.1°C to 37.6°C)	Over 99.6°F (37.6°C)
Rectal	97.3°F to 100.4°F (36.3°C to 38°C)	Over 100.4°F (38°C)
Underarm	95.4°F to 98.6°F (35.2°C to 37°C)	Over 98.6°F (37°C)


*Fever definitions according to the American Medical Association

Changing from Fahrenheit to Celsius

To switch the display between °F and °C and vice-versa: While the thermometer is off, press and hold the On/Scan button  for about 5 seconds until you see the °F or °C. Release button and press it quickly again. Unit has been reset. Wait for display to appear.

Fever Light™

The Fever Light™ occurs during temperature measurement when the temperature is high. (See fever chart included on package or under “Temperature Measurement Range”). If your child has an elevated temperature, red Fever Light™ bar will light up and multiple beeps will sound. High or prolonged fever requires medical attention, especially for young children. Contact your physician.

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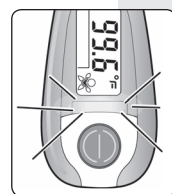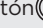
Memory Function

The thermometer memory function stores the last measurement. The last measurement will appear after turning the thermometer on and after it runs its built-in test. This feature allows you to compare your most recent measurement to the next measurement.

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Care and Maintenance

- Clean the base unit and tip with an alcohol swab before and after each measurement.
- Use a soft, dry cloth to clean thermometer body. Never use abrasive cleaning agents, thinners or benzene for cleaning. Do not scratch the surface of the probe lens or the display. Do not expose the thermometer to extreme temperatures, humidity, direct sunlight or shock. Store at room temperature.

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Operating Thermometer

Getting the Thermometer Ready:

- Turn On:** Press and release the On/Scan button . One beep occurs and all segments are displayed with back light for two seconds (Figure 1). The previous measurement is displayed for three seconds (Figure 2). Back light shuts off and “- -” flashes until step 2 is complete.
- Select Mode:** Select oral, rectal, or underarm tip and plug into base unit (Figure 3). The mode is automatically sensed and will display accordingly (Figure 3). **NOTE:** Tip fits in one direction only.
- When you hear beeps, the thermometer is ready to measure temperature.

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Operating Thermometer (continued)

Choosing a Method:

ORAL METHOD—For children 4 years and older.

- Place the probe tip under the tongue in locations indicated by  (Figure 1). Instruct children not to bite down on the probe.

UNDERARM METHOD

- Wipe underarm with a dry towel. Probe must touch skin.
- Place probe tip at center of underarm (Figure 2) and bring arm down against the body.

RECTAL METHOD—Commonly used for babies, young children or when it is difficult to take an oral or underarm temperature.

- Lubricate probe tip with a water-soluble lubricant or petroleum jelly.
- For easier insertion, position baby so his back and legs form a 90-degree (L-shaped) angle. Try placing child on stomach on your lap with his legs down. Gently insert the probe tip no more than 1/2 inch (1.3 cm) into the rectum. Do not force the tip into the rectum if resistance is encountered.

Taking the Temperature:

- Once you have positioned the thermometer, hold it steady while the temperature is obtained. One beep will be heard to indicate that measurement is complete. Back light will light up for 3 seconds. Temperature will be displayed for about one minute without back light. **FEVER ALERT NOTE:** If patient has an elevated temperature, the red Fever Light™ bar will light up and multiple beeps in ten seconds will sound.
- Press the On/Scan button to turn the thermometer off. If left on, automatic shut off will occur in approximately one minute.


NOTE: Wait at least one minute before taking another measurement. Clean the thermometer (see “Care and Maintenance” section).

To Replace Battery

When the “Low Battery” symbol appears (Figure 1), the battery is exhausted and needs replacing.

- Use a coin to loosen battery cover. Remove cover (Figure 2).
- Turn unit upside down and shake battery out.
- Insert a new 3-volt battery (battery number CR1225-3V), correctly positioned with the positive “+” side facing up (Figure 2).
- Place the battery cover on the thermometer and secure the cover.

CAUTION: Replace only with 3V “CR1225” lithium battery. Use of other types of batteries may present a risk of fire or rupture.

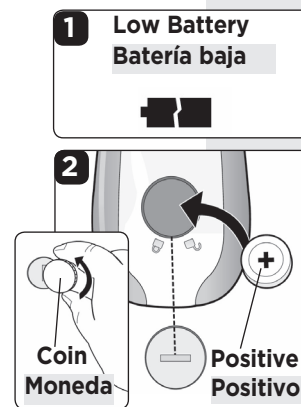WARNING:

- Keep battery out of child’s reach. Swallowing battery could be fatal. Battery should not be charged or placed into extreme heat as it may explode.
- Used batteries should be disposed of properly. Do not dispose of batteries in fire, recharge, put in backwards, disassemble, mix with used or other battery types – they may explode or leak and cause injury.

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©2011 Dorel Juvenile Group. All Rights Reserved. Todos derechos reservados.